

ATTO 647N maleimide

Catalog number: 2857

Unit size: 1 mg

Component	Storage	Amount (Cat No. 2857)
ATTO 647N maleimide	Freeze (< -15 °C), Minimize light exposure	1 mg

OVERVIEW

ATTO 647N is a rhodamine-based fluorescent dye optimized for use in the red spectral region, with similar spectral characteristics as Cy5. It is characterized by high molar absorptivity, a strong fluorescence quantum yield, and excellent thermal and photostability. The dye is moderately hydrophilic and exhibits optimal excitation within the 625-660 nm range, making it compatible with the 647 nm line of Krypton-lon lasers and the 650 nm line of diode lasers. ATTO 647N maintains stable fluorescence across a broad pH range (pH 2-11), supporting its application under diverse experimental conditions. Upon conjugation to a substrate, the dye becomes cationic, carrying a net positive charge of +1. Unlike cyanine dyes, ATTO 647N demonstrates enhanced resistance to atmospheric ozone degradation, which increases its reliability in microarray applications. ATTO 647N is particularly effective for high-precision applications such as single-molecule detection, super-resolution microscopy techniques (e.g., SIM and STED), flow cytometry (FACS), fluorescence in situ hybridization (FISH), and various other biological assays.

The maleimide derivative of ATTO 647N is widely used for labeling biomolecules with free thiol (SH) groups, including antibodies, proteins, thiol-modified oligonucleotides, and low molecular weight ligands. Maleimides react readily with sulfhydryl groups, forming stable thio-ether bonds between the dye and the biomolecule, facilitating robust and reliable labeling for diverse experimental applications.

PREPARATION OF STOCK SOLUTIONS

Unless otherwise noted, all unused stock solutions should be divided into single-use aliquots and stored at -20 °C after preparation. Avoid repeated freeze-thaw cycles

ATTO 647N maleimide Stock Solution (Solution B)

1. Prepare a 10 mM ATTO 647N maleimide stock solution by adding anhydrous DMSO to the vial of ATTO 647N maleimide. Mix well by pipetting or vortexing.

Note: Before starting the conjugation process, prepare the dye stock solution (Solution B) and use it promptly. Prolonged storage of Solution B may reduce its activity. If necessary, Solution B can be stored in the freezer for up to 4 weeks, provided it is protected from light and moisture. Avoid freeze/thaw cycles.

Protein Stock Solution (Solution A)

1. Prepare a 1 mL protein labeling stock solution, by mixing 100 µL of a reaction buffer (e.g., 100 mM MES buffer with a pH ~6.0) with 900 µL of the target protein solution (e.g., an antibody or protein solution with a concentration >2 mg/mL if possible).

Note: The pH of the protein solution (Solution A) should be 6.5 ± 0.5.

Note: Impure antibodies or antibodies stabilized with bovine serum albumin (BSA) or other proteins will not be labeled well.

Note: The conjugation efficiency is significantly reduced if the protein concentration is less than 2 mg/mL. To achieve optimal labeling efficiency, it is recommended to maintain a final protein concentration within the range of 2-10 mg/mL.

Disulfide Reduction (If Necessary)

If your protein does not contain a free cysteine, it must be treated with DTT or TCEP to generate a thiol group. DTT and TCEP are utilized to convert disulfide bonds into two free thiol groups. If using DTT, ensure to remove any free DTT via dialysis or gel filtration before conjugating a dye maleimide to your protein. Below is a sample protocol for generating a free thiol group:

1. To prepare a fresh solution of 1 M DTT, dissolve 15.4 mg of DTT in 100 µL of distilled water.
2. To prepare the IgG solution in 20 mM DTT, first, add 20 µL of DTT stock to each milliliter of the IgG solution while mixing gently. Then, allow the solution to stand at room temperature for 30 minutes without additional mixing. This resting period helps to minimize the reoxidation of cysteines to cystines.
3. Pass the reduced IgG through a filtration column that has been pre-equilibrated with "Exchange Buffer." Collect 0.25 mL fractions as they elute from the column.
4. Determine the protein concentrations and combine the fractions containing the highest amounts of IgG. This can be accomplished using either spectrophotometric or colorimetric methods.
5. Proceed with the conjugation immediately after this step (refer to the Sample Experiment Protocol for details).

Note: IgG solutions should be >4 mg/mL for the best results. The antibody should be concentrated if less than 2 mg/mL. Include an extra 10% for losses on the buffer exchange column.

Note: The reduction can be carried out in almost any buffers from pH 7-7.5, e.g., MES, phosphate, or TRIS buffers.

Note: Steps 3 and 4 can be replaced by dialysis.

SAMPLE EXPERIMENTAL PROTOCOL

This labeling protocol was designed for the conjugation of goat anti-mouse IgG with ATTO 647N maleimide. You may need to further optimize the protocol for your specific proteins.

Note: Each protein requires a specific dye-to-protein ratio, which varies based on the properties of the dyes. Over-labeling a protein can negatively impact its binding affinity while using a low dye-to-protein ratio can result in reduced sensitivity.

Run the Conjugation Reaction

1. Use a 10:1 molar ratio of Solution B (dye)/Solution A (protein) as the starting point. Add 5 µL of the dye stock solution (Solution B, assuming the dye stock solution is 10 mM) to the vial of the protein solution (95 µL of Solution A), and mix thoroughly by shaking. The protein solution has a concentration of ~0.05 mM

assuming the protein concentration is 10 mg/mL and the molecular weight of the protein is ~200KD.

Note: We recommend using a 10:1 molar ratio of Solution B (dye) to Solution A (protein). If this ratio is not suitable, determine the optimal dye/protein ratio by testing 5:1, 15:1, and 20:1 ratios.

2. Continue to rotate or shake the reaction mixture at room temperature for 30-60 minutes.

Purify the Conjugate

The following protocol serves as an example for purifying dye-protein conjugates using a Sephadex G-25 column.

1. Follow the manufacturer's instructions to prepare the Sephadex G-25 column.
2. Load the reaction mixture (From "Run conjugation reaction") to the top of the Sephadex G-25 column.
3. Add PBS (pH 7.2-7.4) as soon as the sample runs just below the top resin surface.
4. Add more PBS (pH 7.2-7.4) to the desired sample to complete the column purification. Then, combine the fractions that contain the desired dye-protein conjugate.

Note: For immediate use, dilute the dye-protein conjugate with staining buffer. If you need to use it multiple times, divide it into aliquots.

Note: For long-term storage, the dye-protein conjugate solution should be either concentrated or freeze-dried.

Characterize the Desired Dye-Protein Conjugate

The Degree of Substitution (DOS) is a key factor in characterizing dye-labeled proteins. Proteins with a lower DOS generally have weaker fluorescence intensity, while those with a higher DOS may also have reduced fluorescence. For most antibodies, the optimal DOS is recommended to be between 2 and 10, depending on the properties of the dye and protein. For effective labeling, the DOS should be controlled to have 5-8 moles of ATTO 647N maleimide per mole of antibody. The following steps outline how to determine the DOS of ATTO 647N maleimide-labeled proteins.

Measure Absorption

To measure the absorption spectrum of a dye-protein conjugate, maintain the sample concentration between 1 and 10 μM . The exact concentration within this range will depend on the dye's extinction coefficient.

Read OD (absorbance) at 280 nm and dye maximum absorption (λ_{max} = 645 nm for ATTO 647N dyes)

For most spectrophotometers, dilute the sample (from the column fractions) with de-ionized water until the OD values fall within the range of 0.1 to 0.9. The optimal absorbance for protein is at 280 nm, while for ATTO 647N maleimide, it is at 645 nm. To ensure accurate readings, make sure the conjugate is free of any non-conjugated dye.

Calculate DOS

You can calculate DOS using our tool by following this link:

<https://www.aatbio.com/tools/degree-of-labeling-calculator>

EXAMPLE DATA ANALYSIS AND FIGURES

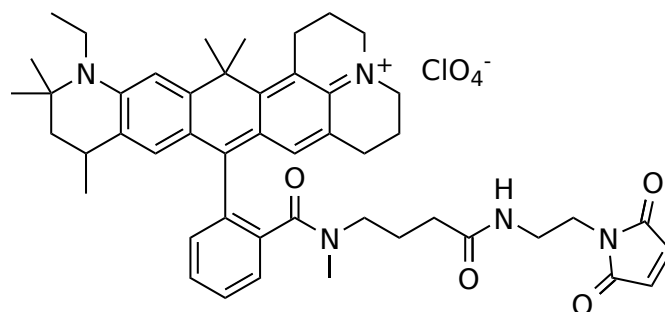


Figure 1. The chemical structure of ATTO 647N maleimide.

DISCLAIMER

AAT Bioquest provides high-quality reagents and materials for research use only. For proper handling of potentially hazardous chemicals, please consult the Safety Data Sheet (SDS) provided for the product. Chemical analysis and/or reverse engineering of any kit or its components is strictly prohibited without written permission from AAT Bioquest. Please call 408-733-1055 or email info@aatbio.com if you have any questions.