

This product is for research use only (not for diagnostic or therapeutic use)

contact: support@agrisera.com

Agrisera AB | Box 57 | SE-91121 Vännäs | Sweden | +46 (0)935 33 000 | www.agrisera.com

Product no AS15 3095

Anti-HEN1 | HUA ENHANCER 1

Product information

Immunogen KLH-conjugated peptide derived from Arabidopsis thaliana HEN1 sequence, Uniprot: Q9C5Q8, TAIR: AT4G20910

Host Rabbit

Clonality Polyclonal

Purity Immunogen affinity purified serum in PBS pH 7.4.

Format Lyophilized

Quantity 50 μg

Reconstitution For reconstitution add 50 μl of sterile water

Store lyophilized/reconstituted at -20°C; once reconstituted make aliquots to avoid repeated freeze-thaw cycles. Please remember to spin the tubes briefly prior to opening them to avoid any losses that might occur from material adhering to

the cap or sides of the tube.

Application information


Recommended dilution 1:1000-1:5000 (WB)

Expected | apparent

104.5 kDa | 105 kDa

Not reactive in No confirmed exceptions from predicted reactivity are currently known

Application example

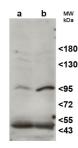

50 μg of total protein from Arabidopsis thaliana whole vegetative rosette wild type Col-0 (a) and HEN1 overexpression mutant (b) extracted with extraction buffer (50 mM Tris pH 7.5; 150 mM NaCl; 1 mM EDTA; 10% v/v Glycerin; 1 mM DTT, 1x Complete Protease Inhibitor Cocktail, Roche) and denatured with Laemmli buffer at 95 °C/5 min., were separated on 7.5% SDS-PAGE and blotted 1.5 h to PVDF using tank transfer. Blots were blocked with blocking buffer (5% milk powder; 1x TBS; 0.1% Tween-20) overnight at 4°C with agitation. Blot was incubated in the primary antibody at a dilution of 1:1000 for 2h at RT with agitation. The antibody solution was decanted and the blot was rinsed briefly and then washed tree times for 15 min in TBS-T at RT with agitation. Blot was incubated in secondary antibody (anti-rabbit IgG horse radish peroxidase conjugated, from Agrisera, AS09 602) diluted to 1:20 000 in blocking buffer for 1h at RT with agitation. The blot was washed as above and developed for 5 min with Amersham ECL Prime and expose to Amersham Hyperfilms ECL for 20 seconds.

Courtesy of Dr. Pablo Manavella, Instituto de Agrobiotecnología del Litoral (IAL), Argentina