

Separase rabbit pAb

Cat No.:ES3423

For research use only

Overview

Product Name	Separase rabbit pAb
Host species	Rabbit
Applications	WB;IHC;IF;ELISA
Species Cross-Reactivity	Human;Mouse
Recommended dilutions	Western Blot: 1/500 - 1/2000.
	Immunohistochemistry: 1/100 - 1/300.
	Immunofluorescence: 1/200 - 1/1000. ELISA:
	1/40000. Not yet tested in other applications.
Immunogen	The antiserum was produced against synthesized
	peptide derived from human SEPARASE. AA
	range:767-816
Specificity	Separase Polyclonal Antibody detects endogenous
	levels of Separase protein.
Formulation	Liquid in PBS containing 50% glycerol, 0.5% BSA and
	0.02% sodium azide.
Storage	Store at -20°C. Avoid repeated freeze-thaw cycles.
Protein Name	Separin
Gene Name	ESPL1
Cellular localization	Cytoplasm. Nucleus.
Purification	The antibody was affinity-purified from rabbit
	antiserum by affinity-chromatography using
	epitope-specific immunogen.
Clonality	Polyclonal
Concentration	1 mg/ml
Observed band	230kD
Human Gene ID	9700
Human Swiss-Prot Number	Q14674
Alternative Names	ESPL1; ESP1; KIAA0165; Separin; Caspase-like
	protein ESPL1; Extra spindle poles-like 1 protein;
	Separase
Background	Stable cohesion between sister chromatids before
	anaphase and their timely separation during
	anaphase are critical for chromosome inheritance.
-	



+86-27-59760950


ELKbio@ELKbiotech.com

www.elkbiotech.com

23-2, No.388 Gaoxin 2nd Road, Wuhan East Lake Hi-tech Development Zone, Hubei , P.R.C




In vertebrates, sister chromatid cohesion is released in 2 steps via distinct mechanisms. The first step involves phosphorylation of STAG1 (MIM 604358) or STAG2 (MIM 300826) in the cohesin complex. The second step involves cleavage of the cohesin subunit SCC1 (RAD21; MIM 606462) by ESPL1, or separase, which initiates the final separation of sister chromatids (Sun et al., 2009 [PubMed 19345191]).[supplied by OMIM, Nov 2010],

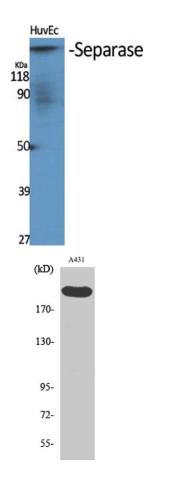

Western Blot analysis of various cells using Separase Polyclonal Antibody diluted at 1:1000

Western Blot analysis of A431 cells using Separase Polyclonal Antibody diluted at 1:1000

Immunofluorescence analysis of HUVEC cells, using SEPARASE Antibody. The picture on the right is blocked with the synthesized peptide.

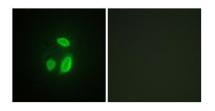



+86-27-59760950


ELKbio@ELKbiotech.com

www.elkbiotech.com

23-2, No.388 Gaoxin 2nd Road, Wuhan East Lake Hi-tech Development Zone, Hubei , P.R.C



Immunohistochemistry analysis of paraffin-embedded human colon carcinoma tissue, using SEPARASE Antibody. The picture on the right is blocked with the synthesized peptide.

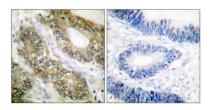



+86-27-59760950

ELKbio@ELKbiotech.com

www.elkbiotech.com

23-2, No.388 Gaoxin 2nd Road, Wuhan East Lake Hi-tech Development Zone, Hubei , P.R.C